



# Double Cell Coil User Manual






These instructions refer to the use of the coils with a cell phone but any mobile device (tablet, laptop, etc.) may be used.

## Getting Started

1. Remove the coils and cable from the box and remove the plastic bag on the cable.
2. If your cell phone does not have a circular audio jack, an audio adaptor must be purchased and connected to the cell phone.
3. Insert the cable connector into the audio jack on the cell phone or the audio adaptor.






4. Scan this QR code with the cell phone camera to go the YouTube channel playlist or follow steps 5 and 6.
5. Open your web browser and type in cellcoil.com and touch the MENU icon.
6. Touch the THERAPY PROGRAMS option.
7. Choose the program you want and press the PLAY button.

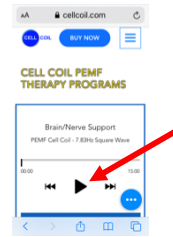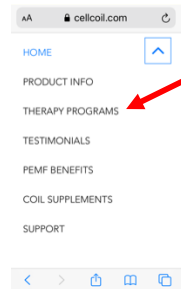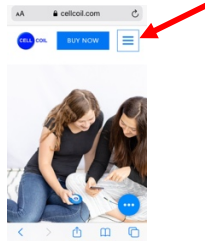

8. Turn the volume up to 100%.
9. To pause or stop the therapy press the PAUSE button.
10. To start the therapy again press the PLAY button.
11. The therapy will continue to play even when the screen is off.

The audio therapy will start to play and the coil will radiate a magnetic field at the specified frequency.

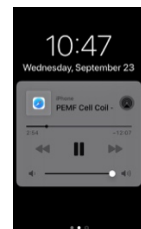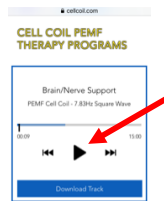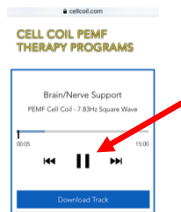

12. Place the coils on the areas to be treated and hold them until the therapy is done, or use the covers and strap to hold the coils in place (see the following steps).
13. Remove the covers and the strap from the product box.

**Always place the Cell Coil logo sides against the areas to be treated.**



14. Fold the cover flaps back and insert a coil in each cover.

Make sure to insert the coils so that the Cell Coil logo sides will be against the areas to be treated.



15. Fold the cover flaps over the coils and attach them to the small Velcro strips.



16. Unwind the strap. The strap may be cut to a convenient length and reconnected a later time.



17. Place the coils on the areas to be treated with the Velcro strips on the covers to hold them in place until the therapy is done.





18. Attach the strap to the Velcro strips on the covers to hold them in place until the therapy is done.



## When The Therapy Is Finished

1. Stop the therapy by touching the Cell Coil icon and pressing the PAUSE button.
2. Remove the coils and strap from the area being treated.
3. Remove the coil cable connector from your phone.

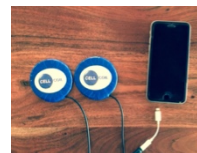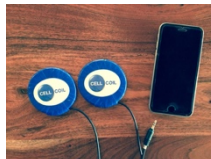

## Support

Scan this QR code with the cell phone camera to view a video about how to use the cell coils.

Website: [cellcoil.com](http://cellcoil.com)

Support: (844) 448-1616



Cell Coil products have not been evaluated by the FDA. The Cell Coil is intended for the beneficial effects as set forth in the directions and instruction literature. This product does not claim to diagnose, treat, cure, or prevent any medical condition. This device should not be used in the following conditions: pregnancy, epilepsy, bleeding, implanted electronic devices or metal implants. Use as directed by your health care professional.